

Operation Manual

Shodex MSpak GF-310 4D

(Please read this operation manual carefully to achieve the best and consistent column performance for a long time.)

Important Handling Instructions

Caution!

- Please consult the Safety Data Sheet (SDS) of reagents and solvents used with the column and understand their proper handling methods to prevent potential health hazards or death from occurring.
- Please wear appropriate personal protective equipment such as lab goggles and gloves when handling organic solvents and acid and alkaline reagents. Avoid any direct physical contact to prevent chemical injuries.

Before Using the Column

- (1) Please visually inspect the column package and the column surface for any damage.
- (2) Please check the product name and serial number (Serial no. or S/N) written on the column package and adhesive label on the column body.
- (3) Please download the Certificate of Analysis (CoA) for the purchased product. The CoA can be downloaded from Shodex website (<https://www.shodex.com/download/>). You will be asked to enter the serial number.

1. Introduction

Thank you for purchasing the Shodex product. Shodex MSpak GF-310 4D is a polymer-based size exclusion (SEC) chromatography columns with high solvent tolerance. In addition to the SEC analysis, it can be used for LC/MS analysis of drugs in biological samples.

2. Column Components

Please refer to the Shodex website: <https://www.shodex.com/en/da/07.html>

3. Column Specifications

Product Code	Product Name	Column Size (mm)		Particle Size (µm)	Theoretical Plate Number (Per Column)	Target Molecular Weight Range*	Exclusion Limit*
		I.D.	Length				
F7600110	MSpak GF-310 4D	4.6	150	5	≥ 10,000	300 - 30,000	40,000

*1Reference value only / Measured with pullulan

Base Material : Spherical porous particles of polyvinyl alcohol

Column Housing : SUS-316

Screw Type : Internally-threaded type No.10-32 UNF

Shipping Solvent : Water

4. Usable Conditions

4.1 System Settings

Product Name	Flow Rate (mL/min)		Maximum Pressure (MPa/column)	pH Range	Temperature Range (°C)
	Recommended	Maximum			
MSpak GF-310 4D	0.2	0.5	4.0	2 - 9	4 - 60

Usable solvents are listed below.

- (1) The standard eluent is a mixture of water and acetonitrile or water and methanol of any ratio.
- (2) Buffers and aqueous solutions of different salts can be used separately or together. Usable buffers include phosphate, acetate, citrate, and Tris. Usable aqueous salt solutions include sodium chloride, potassium chloride, sodium sulfate, potassium sulfate, and ammonium sulfate. Please keep total concentration of salts under 0.1M.

Attention!

- Use the column within above stated flow rate, pressure, and temperature ranges. Using the column outside the given range may damage the column and lower its performance.
- When using a mixture of buffer (or aqueous solution of salt) and organic solvent, make sure there is no precipitation of salt.
- When using highly corrosive salts such as sodium chloride, wash out the salts at the end of analysis. The metal parts of the devices and/or the columns may rust.
- Column pressure is influenced by eluent composition, flow rate, and column temperature. When changing the eluent compositions, adjust the flow rate and column temperature so that the column pressure remains below the usable maximum pressure.

5. Eluent Preparation

- (1) Degas the eluent fully to prevent the formation of air bubbles.
- (2) Presence of small debris or insoluble substances may result in deterioration of columns and/or they may appear as noise on chromatograms. Filter the eluent with a 0.45- μ m disposable filter to prevent the problems from occurring.

Attention!

- Whenever water is required, use ultra-pure water freshly generated by a water purification system or water from a newly opened HPLC grade distilled water bottle. Use of HPLC grade organic solvents of guaranteed quality, which can be used without problems in HPLC is recommended. If organic solvents with different grades are used together, make sure that their qualities are all suitable for the analysis prior to the use. Solvents left in opened bottles for a long time should not be used. The content may have been changed, absorbed moisture, or has been contaminated.
- Always use freshly prepared solvents. Solvents stored for a long time may have changed their compositions and may influence elution patterns and/or damage the column.

Note

- Use of an on-line degasser is recommended.

6. Sample Preparation

- (1) If possible, use the eluent for analysis to dissolve or dilute samples. If this is difficult, use a solvent which has a composition that is as close as possible to the eluent composition and which fully dissolves or dilutes the sample.
- (2) Filter diluted sample solutions using disposable 0.45- μ m filters to prevent the column from clogging or deteriorating.
- (3) Suggested injection volume is less than 20 μ L.

Attention!

- When a sample is dissolved in a solvent other than the eluent and if the sample matrix contains components which do not dissolve in the eluent fully, precipitates may form and clog the column.

7. Column Usage Procedure

MSpak GF-310 4D is a semi-micro column. To achieve its best performance, please use it with a semi-micro type HPLC system.

7.1 HPLC System Preparation


Wash the entire HPLC system prior to column installation, including all flow-lines and sample loop by switching the valves, and then replace the washing solution with the eluent to be used. If desired new eluent has low miscibility/solubility to the eluent of previous analysis, first use the eluent that is miscible/soluble to both eluents, and then replace it with the desired eluent.

Attention!

- If the eluent left in the HPLC system is not compatible with the column to be used, it may damage the column.
- A drastic change in the eluent compositions may remove substances adsorbed on the HPLC system and they may enter and deteriorate the column.

7.2 Column Installation

- (1) Connect the column to HPLC system by following the “flow direction arrow” (→) indicated on the column adhesive label. If a guard column is used, position the guard column in front (before the inlet) of the analytical column.
- (2) Make sure to insert the tubing all the way to the end fitting and secure it with the male nut. It is important that there is no extra space between the tubing and the column side of the end fitting. Presence of an extra space will let the sample to spread out and may result in wide peaks.

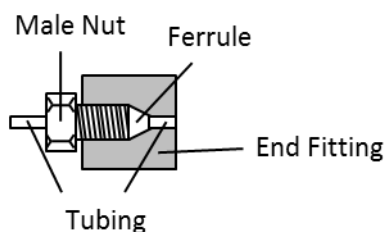

- (3) Set the initial flow rate to less than 0.1 mL/min and start the system. If the column is to be heated during the analysis, keep the low flow rate until the column temperature reaches to the set temperature, and then gradually increase the flow rate to the desired temperature.

Caution!

- Verify that there is no solvent leak. The solvent leak may cause electronic leakage, rust, and/or chemical injury.

Attention!

- Make sure not to let air bubbles enter the column while installing the column. The air bubbles may damage the column.
- When restarting the system after column installation or after holding the eluent flow, start the system at less than 0.1 mL/min flow rate. A rapid increase in pressure can damage the column.
- If the column was heated during the analysis, lower the flow rate to less than 0.1 mL/min at the end of analysis. Then, turn off the column oven to let the column temperature returns to room temperature before stopping the pump. This is to prevent creating an empty space in the column, which deteriorates the column. Since if the pump was stopped while the eluent inside the column is still hot, the eluent volume decreases and creates an empty space when the eluent temperature decreases.

Note

- It is recommended to set the pump limiter to avoid exceeding the maximum pressure.

7.3 Solvent Exchange

To replace the solvent, set the system at less than 0.15 mL/min. Recommended solvent volume to introduce at each step is 3 to 5 times of the column volume.

- (1) Check miscibility/solubility of the desired new solvent and the solvent currently filled in the column.
- (2) When replacing the current solvent with a solvent with low miscibility/solubility to the current solvent, first use a solvent that is miscible/soluble to both solvents, and then replace it with the new solvent.
e.g. When replacing from a highly concentrated buffer solution to water/acetonitrile, first run water and then replace it with the eluent.
- (3) When using a gradient method, changes in the eluent compositions may increase the column backpressure. Adjust the flow rate and column temperature so that the column backpressure remains below the usable maximum pressure throughout the analysis.

7.4 Column Cleaning

Problems in peak shapes and elution time changes or elevated column pressure are often caused by the deposition of insoluble or adsorbing components from the sample/flow-line inside the column. These problems may be resolved by cleaning the column.

In case multiple number of analytical columns are used together, wash them separately. During the column cleaning, disconnect the detector and collect the washing solution directly from the column outlet into a waste container (i.e., do not let the solution go through the detector).

If the column performance does not improve (recover) after performing the column cleaning, please replace the column with a new one.

<Cleaning method>

- (1) Insoluble components that block the column inlet may be removed by reversing the flow direction, i.e., introducing the eluent from the column outlet, with flow rate at less than half of the recommended flow rate.
- (2) Follow below cleaning steps for adsorbing components. For an efficient cleaning, reverse the flow direction. Set the flow rate at less than 0.15 mL/min. Recommended solvent volume to introduce is 5 to 10 times of the column volume.

Method 1: Adsorption of hydrophobic compounds

Introduce the eluent with higher acetonitrile or methanol concentration.

Method 2: Adsorption of ionic compounds

Introduce the eluent with a higher salt concentration.

Attention!

- Keep the salt concentration within the concentrations stated in 4. Usable Conditions.
- Keeping the washing solution in the column for a long time will lead to column deterioration. Please replace the washing solution with the eluent immediately after cleaning.

8. Column Storage

Remove the column from HPLC system after replacing the in-column solvent with the initial shipping solvent. Securely tighten the end caps and store the column at a location with stable temperature (a cool and dark space is recommended). Refer to section 7.3 Solvent Exchange for how to replace the eluent.

Attention!

- Never allow inside the column to dry. It can damage the column.

9. Column Inspection

Please refer to the inspection method described in the CoA. At Shodex, “half width method” is adopted for the calculation of plate count and “asymmetry factor (Fas)” is adopted for the calculation of peak symmetry. Please refer to the Shodex website for the detail: <https://www.shodex.com/en/da/07.html>

Attention!

- Plate count and Fas values change significantly depend on samples and/or analysis conditions being used. To check the initial column condition, please make sure to use the same sample and the analysis condition mentioned in the CoA.

10. Additional Warnings

- (1) Do not remove end fittings.
- (2) Do not make a strong impact on the column. Do not drop or hit the column on a hard surface.
- (3) Please follow a proper waste disposal method specified by your local regulations.

Please refer to the Shodex website (<https://www.shodex.com/>) for product details and their applications. For additional assistance, contact the distributor from whom you purchased the column or contact your regional Shodex support office (https://www.shodex.com/en/support_office/list).